

BCI practical course : “Hello World” & ERP Viewer

Jason Farquhar

<J.Farquhar@donders.ru.nl>

Learning Goals

Understand:

- What is needed to make a BCI, i.e. progress tracking, data acquisition, annotation and processing, stimulus presentation, and an overall process scheduler/sequencer
- How to use event-driven programming ideas coupled to a global shared event pool (blackboard) to provide these facilities
- How the fieldtrip buffer provides the event blackboard which is used for inter-process communication.

Know how to:

- What the struct of an 'event' is and how to use it to annotate data with experiment relevant event information
- present simple visual stimulus/feedback to the user/experimenter
- How to wait for specific events, get the necessary data, process it and post the updated results back to the event blackboard
- Startup the buffer and an experiment control Matlab, and how to connect these processes to provide a basic BCI
- Test your experiment with simulated data generated by the signal-proxy
- Debug your experiment when it fails!

Today's Plan

- Discussion : What do we need to make a BCI?
- Introduction to the Buffer-BCI framework

break

- Hands-on 1: Echo-server (event IPC)
- Hands-on 2: Hello World
- Hands-on 3: Sequenced Sentences


break

- Hands-on 4: Visual ERP Viewer


Discussion: What do we need to make a BCI?

- . Based on your prior knowledge and experience with the hands on demo we've did last time.
- . Discuss: What do we require to make a BCI system?
- . Think about:
 - Hardware requirements?
 - Software requirements?
 - Information flows?

BCI information flow

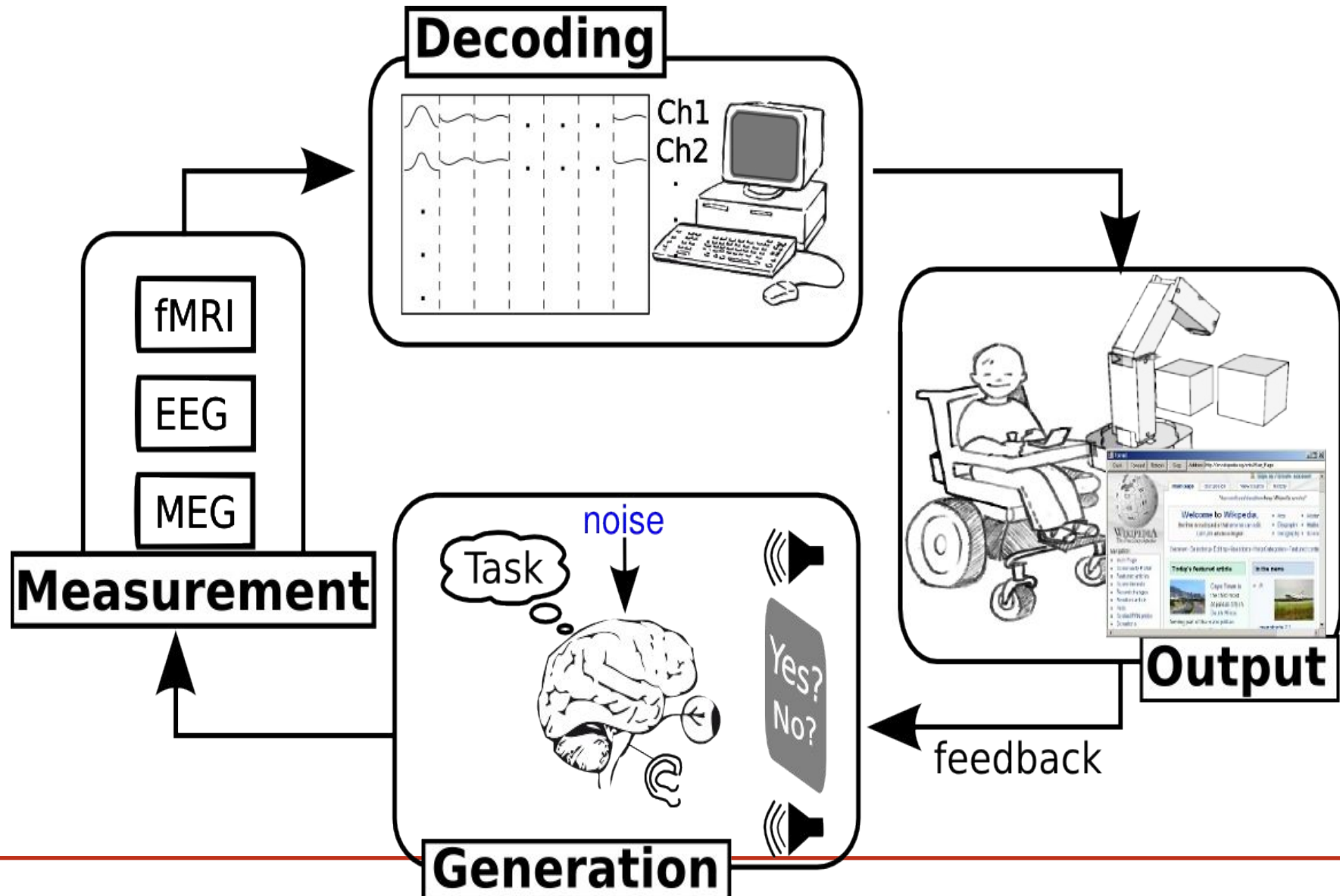


Gross structure of a typical BCI experiment

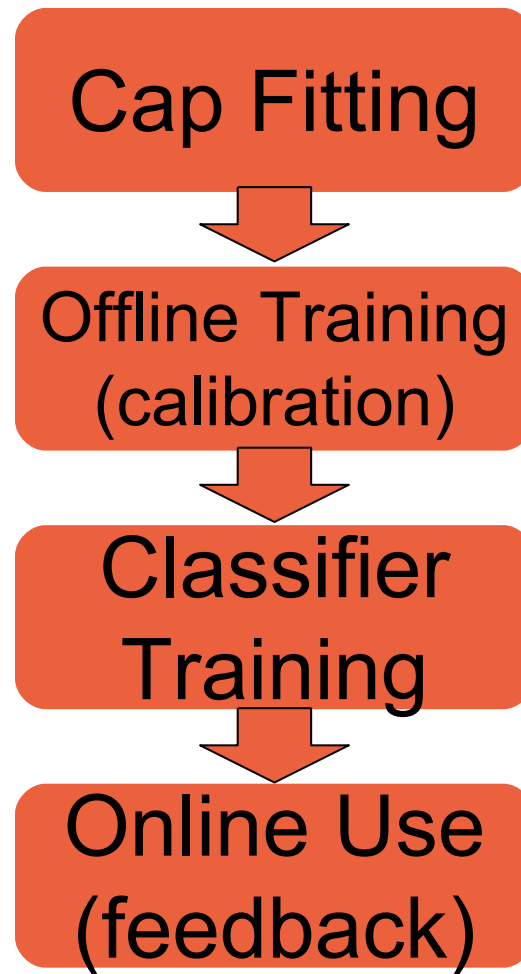


BCI terminology (our group!)

- . Epoch/Trial
 - . Single BCI prediction
 - . e.g. 1-imagined movement, 1 visual-speller flash
- . Sequence
 - . Short group of epochs (~1min)
 - . v. short breaks 1-2sec between epochs (usually automatic)
 - . short (usually self-paced) subject break between sequences (~10sec)
- . Block/Run/Phase
 - . Short group of sequences (>10min)
 - . long (~1-2min) subject break between blocks
 - . e.g. cap-fitting, calibration, classifier training, on-line use
- . Session – during one cap-fitting
- . Experiment – imagined movement, visual-speller etc.


Flow chart of an individual epoch in a simple BCI experiment

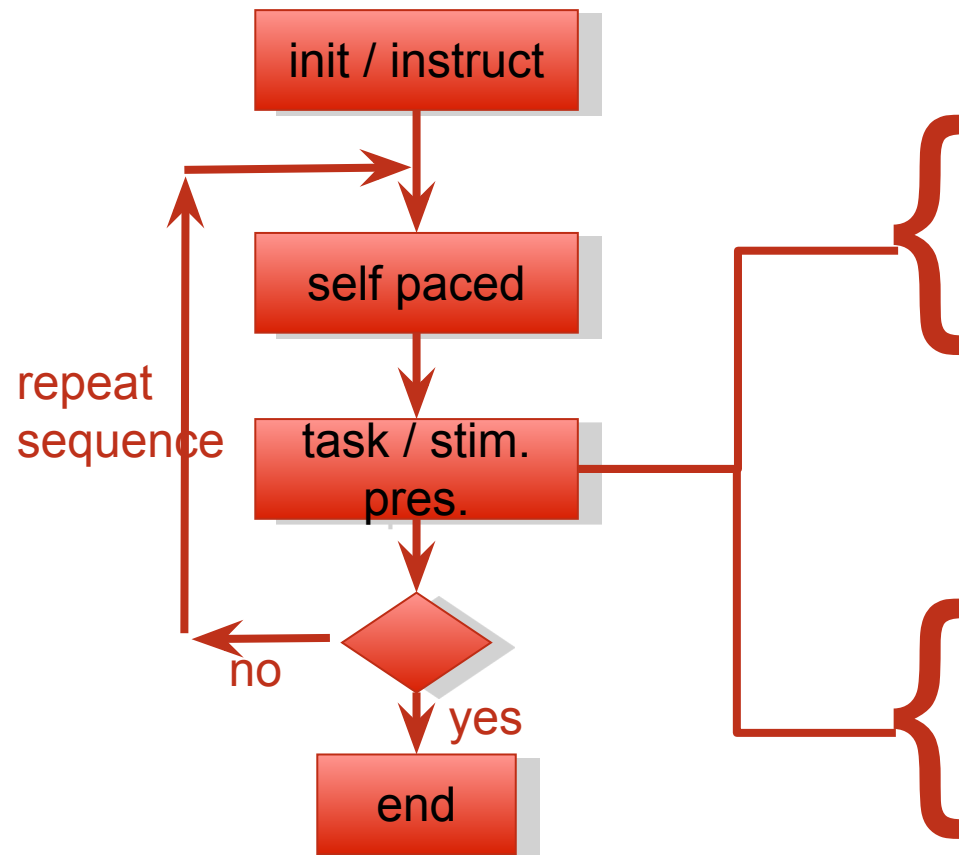

Processing epochs offline

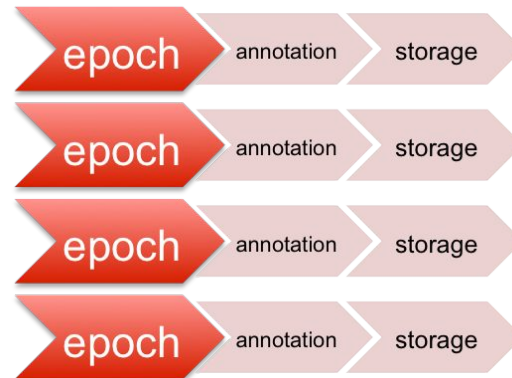

Processing epochs online

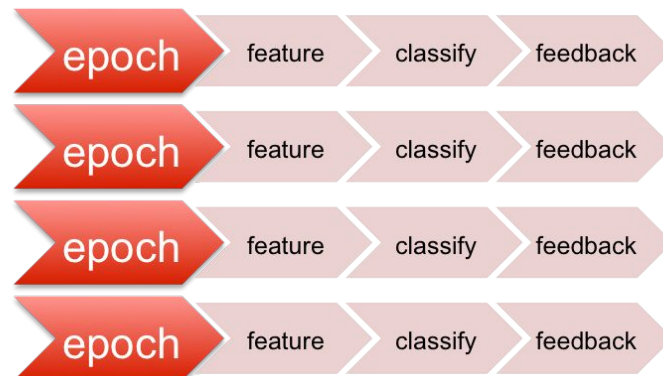


Requirements: what do you need to build a BCI?

- 1) Way of **tracking** where we are in execution of the experiment flowchart, i.e. block, sequence, epoch number.
- 2) Way of **annotating** data to what the subject was experiencing/doing at that time with what was measured from their brain/body, e.g. LH movement, reading instruction, watching queue, etc.
- 3) **Data acquisition**: Drivers to extract data from hardware (and combine data from different hardware sources)
- 4) **Stimulus Generation**: makes stimuli that the subject will experience, for subject instruction, feedback, event-related stimuli
- 5) Something to **process** the signals, firstly to train the classifier, and secondly to decode the users mental state, i.e. do the BCI bit ;-)
- 6) **Scheduler** (sequencer?) to tie it all these bits together,
 - so the correct functions, i.e. stimulus display, signal processing, are executed
 - at the correct position in the experiment flowchart
 - based on the right bits of measured data

Summary


- To build a BCI we need a system to; **track** our progress through the experiment, **acquire**, **annotate** and **process** data, present the **stimuli** and **schedule** all these processes in an appropriate way.
- Next we introduce a Matlab based system which provides these facilities.

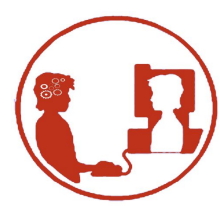

Overview – “buffer_bci” framework

- Open-source BCI development environment
- Available from Github:
www.github.com/jadref/buffer_bci
- Core is “[fieldtrip-buffer](#)”
 - Hardware access, data-storage, IPC
 - Platform independent: Win, Mac, Linux, (Android/iOS)
 - Language Neutral: C, C++, Matlab, Python, C#, java
- Matlab/Octave/Python/Java based BCI examples



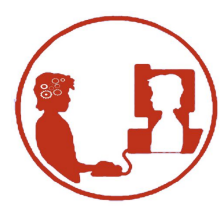

Overview – “buffer_bci” framework

- Open-source BCI development environment
- Available from Github: www.github.com/jadref/buffer_bci
- Demos:
 - Brain Viewer
 - Games : Snake, Pacman, Sokoban
 - Visual Speller
 - Thought-buttons
- Tutorial:
 - How a BCI works and how to build one



Buffer-BCI Framework

.We can break the requirements into 4 largely independent **communicating** processes:

1)Data-acquisition & annotation

- .Get data from hardware
- .Attach annotations (markers, events) to particular data sample

2)Experiment control (scheduling)

- . Control the flow of the experiment

3)Stimulus generation

- . Make stimuli when requested by the expt controller
- . Make feedback based on predictions generated by the sig-processor


4)Signal Processing

- .Process the data based on the annotations, and generate predictions

Buffer-BCI Framework

Basic Idea:

- set of independent processes
- any process can send/recieve **data-annotation events**
- events are visible to **all** other processes
- Processes **communication** implemented by sending recieving events
- (N.B. As all events are saved with the data, annotations are automatically archived for later off-line use).
- Similar in concept to that used in 'Blackboard architectures' for AI, see en.wikipedia.org/wiki/Blackboard_system

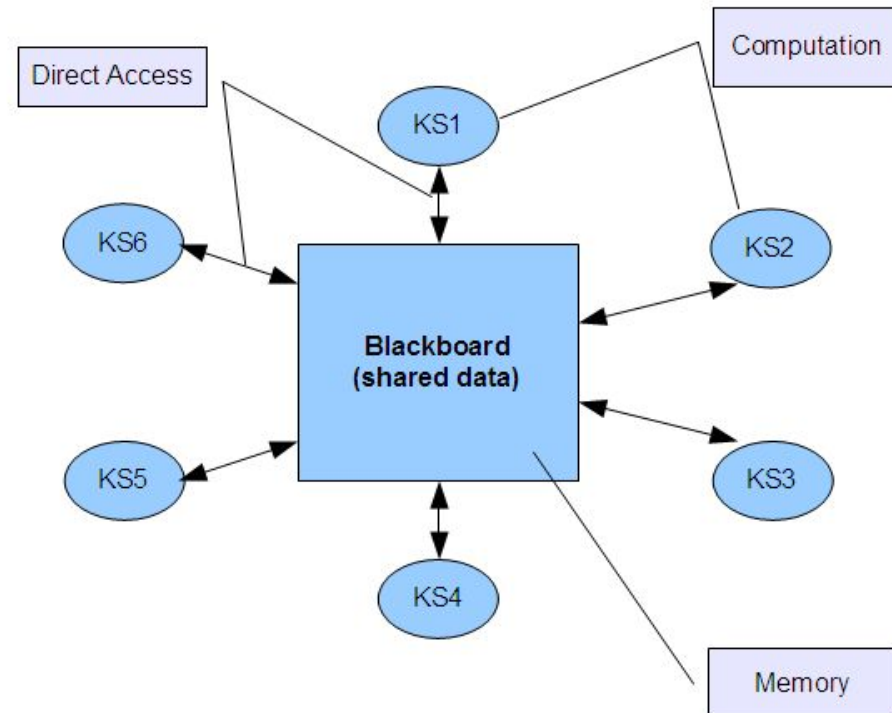


Ft-buffer based Implementation

- Buffer-BCI framework implemented using the fieldtrip-buffer system (fieldtrip.fcdonders.nl/development/realtime)
- Ft-buffer provides:
 - Drivers for data-acquisition
 - 1) buffer storage for **data** (~last 1 minute data)
 - 2) buffer storage for **events** (~last 50 events)

Idea:

- Buffer events store represents the **blackboard** used for inter-process communication (IPC)
- Every event has **timestamp** (sample number) used for data-annotation

Typical Structure of BCI system

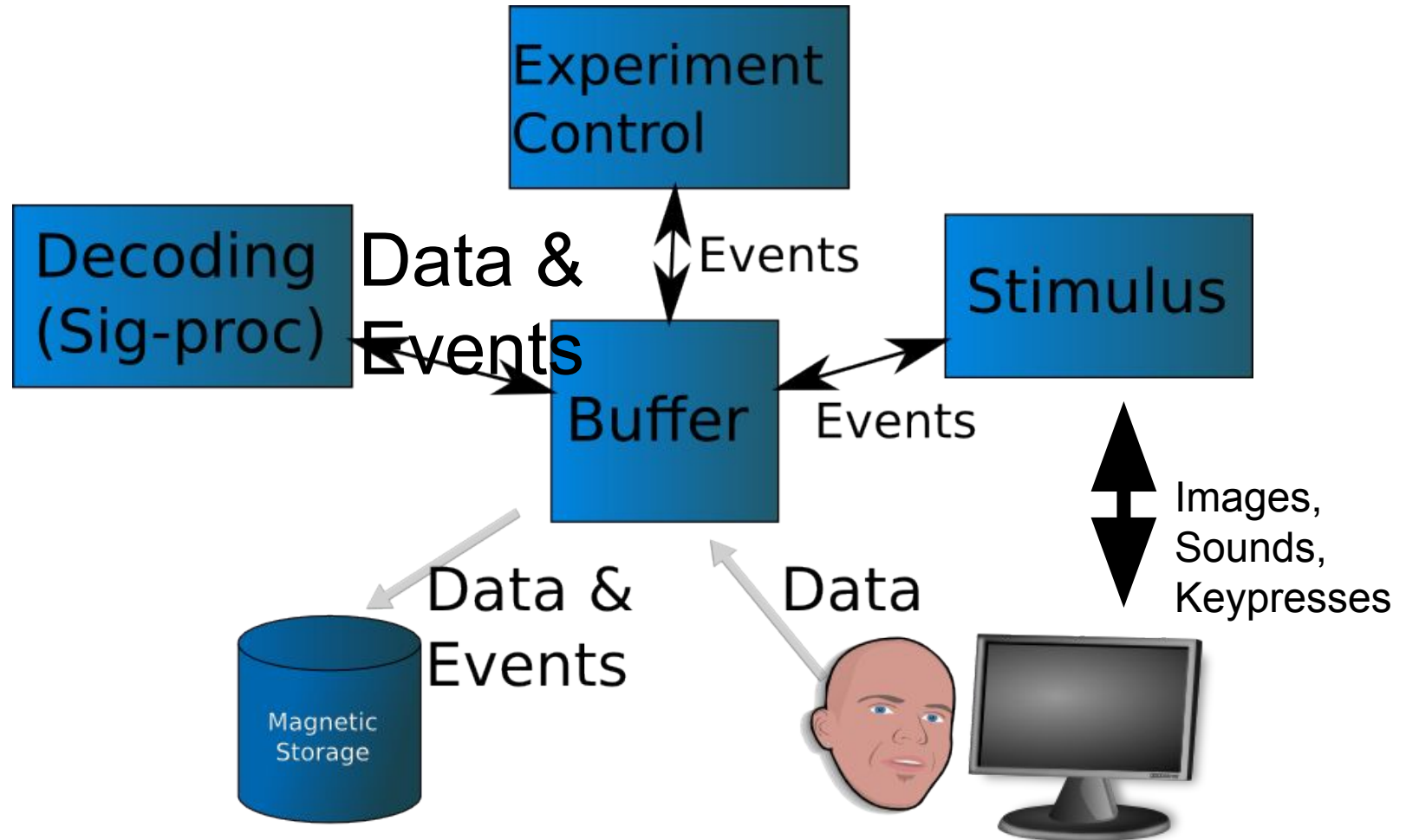






Buffer_bci Usage

client/server architecture

- 1) **Buffer server:** for data storage, annotation (IPC) and archive
- 2) **Acquisition driver:** to access (or simulate) the measurement hardware
- 3) **BCI Application.**
- 4) Normally, split into 2/3 pieces:
 - 1) **Signal-processing**
 - 2) **Stimulus** presentation
 - 3) **Experiment control**
 - 4) (commonly part of stimulus-presentation)

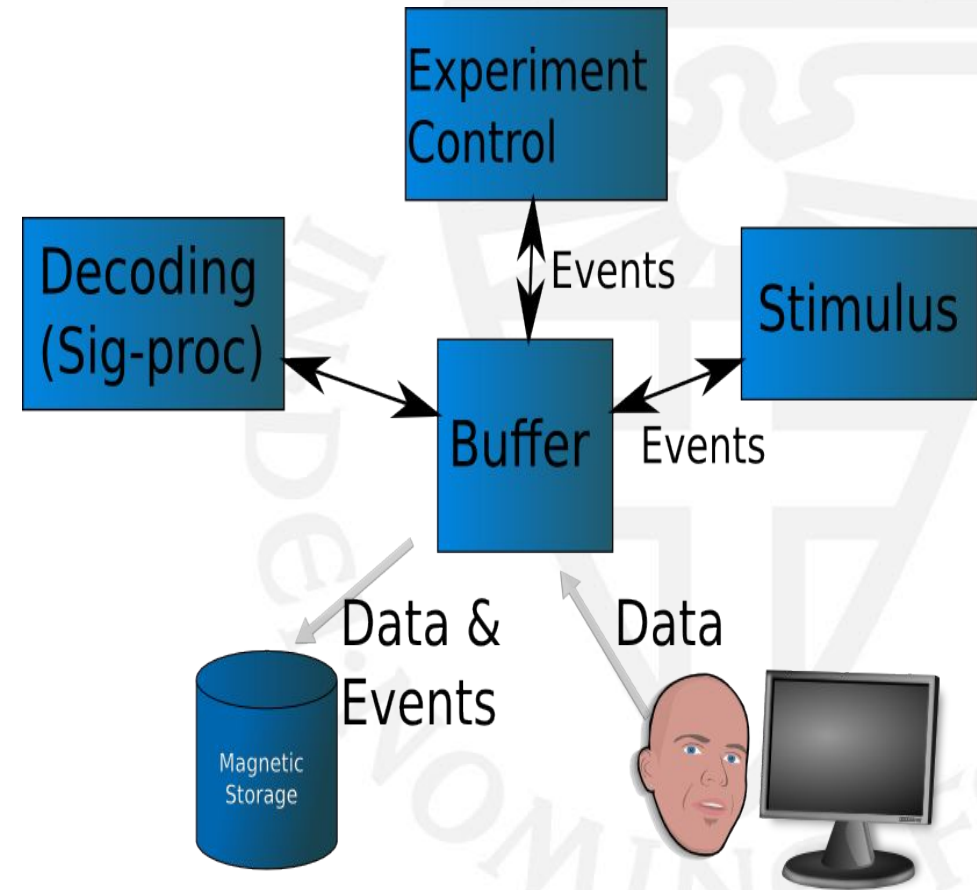

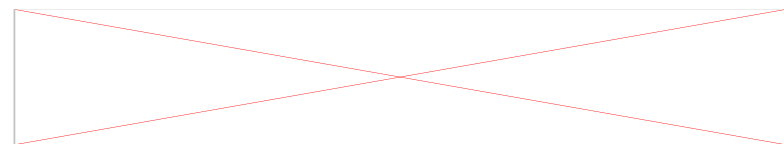


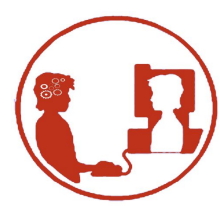

Buffer_bci Quickstart

client/server architecture

- 1) **Buffer server:** for data storage, annotation (IPC) and archive
- 2) **Acquisition driver:** to access (or simulate) the measurement hardware
- 3) **BCI Application.**
- 4) Normally, split into 2/3 pieces:
 - 1) **Signal-processing**
 - 2) **Stimulus** presentation
 - 3) **Experiment control**
 - 4) (commonly part of stimulus-presentation)

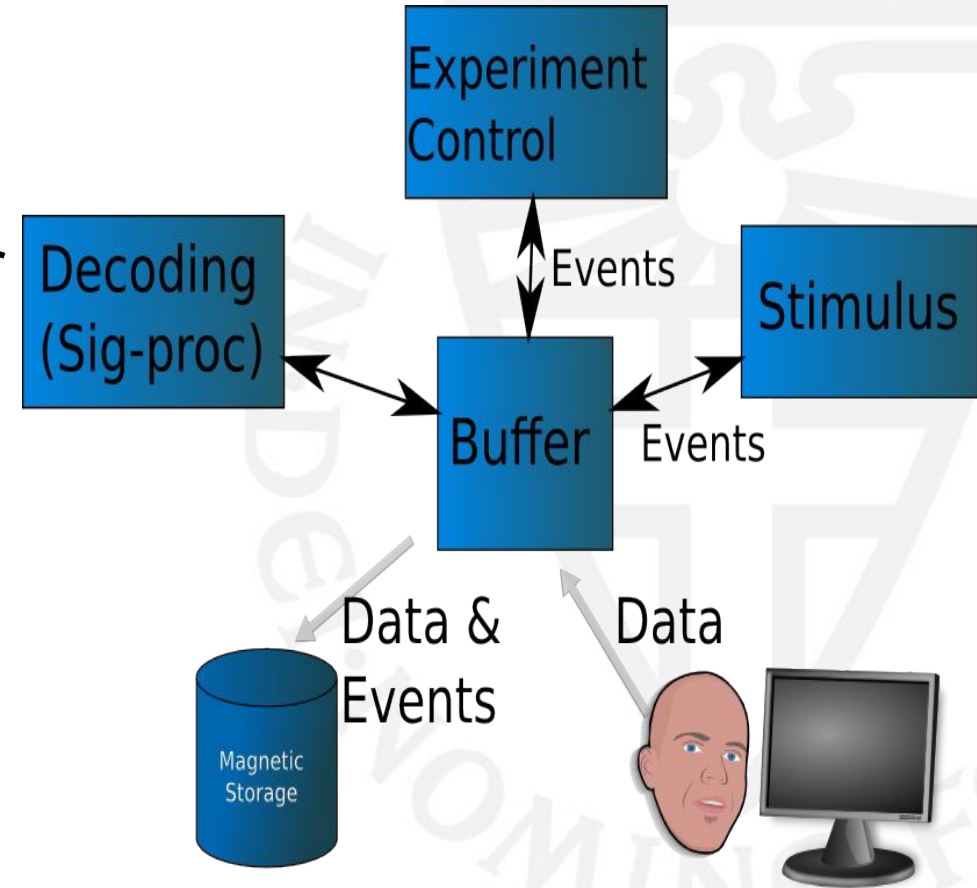

Running buffer-bci code:

You need to have (**at least**) the following processes running:

- 1) Buffer server
- 2) Data-acquisition system (real or simulated)
- 3) BCI Application

1) Buffer Server

- Saving server (saves all data and events)
- `dataAcq/startJavaBuffer.{sh/bat}`
- Server (saves nothing):
`dataAcq/startJavaNoSaveBuffer.{sh/bat}`

C-based implementations:

- Saving: `dataAcq/startBuffer.{sh/bat}`
- NoSaving: `dataAcq/startNoSaveBuffer.{sh/bat}`

2) Data-acquisition

Simulated data: (for debugging/testing)

- **Java-based:** `dataAcq/startJavaSignalProxy.{sh/bat}`
- (Alternative implementations Matlab: `startMatlabSignalProxy`, c: `startSignalProxy`)

Real data: (for debugging/testing)

- **FilePlayback:** `dataAcq/startJavaFileProxy.{sh/bat}`
- **Microphone:** `dataAcq/startAudio.{sh/bat}`

EEG amplifier data: (for real usage)

- **Biosemi:** `dataAcq/startBiosemi.{sh/bat}`
- **Mobita:** `dataAcq/startMobita.{sh/bat}`
- **OpenBCI:** `dataAcq/startOpenBCI.{sh/bat}`
- **Many others...** (look in `dataAcq/startXXXXX.{bat/sh}`)

3) BCI-Application

- Generally what you will write yourself in this course ;-)
- Typically in 2 (or more) separate processes:

1) Stimulus presentation:

- Controls what the user sees/hears, i.e. UI
 - Show App, show feedback, control output-devices, etc.

2) Signal processing:

- Processes the EEG data and events to generate predictions for UI feedback.
 - Collect calibration data, train classifier, generate predictions

Example : Games-demo

1) Buffer Server:

1) `dataAcq/startJavaNoSaveBuffer.{sh,bat}`

2) Data Acquisition (simulated) :

1) `dataAcq/startJavaSignalProxy.{sh,bat}`

3) BCI Application(1) -- Stimulus Presentation:

1) `games/runGame.{bat,sh}`

4) BCI Application(2) – Signal Processing:

1) `games/startSigProcBuffer.{bat,sh}`

Quickstart Example : Games-demo

1) Core Quickstart: Server + Acquisition + EventViewer:

1) Debug mode (simulated data):

2) `debug_quickstart.{sh,bat}`

3) EEG Mode (mobita source):

4) `eeg_quickstart.{sh,bat}`

2) BCI Application(1) -- Stimulus Presentation:

1) `games/runGame.{bat,sh}`

3) BCI Application(2) – Signal Processing:

1) `games/startSigProcBuffer.{bat,sh}`

Useful (debugging) Functions: EventViewer

Seeing the what events are sent and when is important for debugging experiment flow.

Basic event viewer has been implemented in multiple languages

MATLAB:

- **Quickstart:** `dataAcq/startMatlabEventViewer.{sh,bat}`
- **Source code:** `matlab/utilities/eventViewer.m`

JAVA :


- **Quickstart:** `dataAcq/startJavaEventViewer.{sh,bat}`
- **Source code:** `java/echoClient/eventViewer.java`

Python:

- **Source code:** `python/echoClient/eventViewer.py`

C :

- **Source code:** `c/echoClient/eventViewer.c`

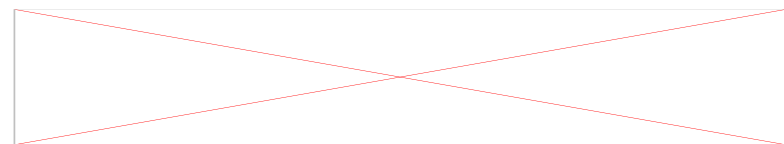

Hands-on 1: Event Echo-Server

Events for IPC

- As well as being used for data annotation, Events are used for inter-process-communication,
 - e.g. to communicate results from signal-processing to stimulus presentation
- To use events in this way, each process must
 - **monitor** for new events
 - **filter** out the events it should react to
 - **send** response events

(Key concept) event structure

sample

- . time at which event occurred in samples from start of experiment

type

- . arbitrary event type (usually a string)

value

- . arbitrary event value (usually string or number)

duration (optional)

- . duration of the event in samples

offset (optional)

- . zero-time for the event.
- . Usually, offset from sample at which the event actually started.

Examples:

Visual speller “flash”;

```
ev=struct('sample',123,...  
         'type','stimulus.flash',...  
         'value',[0 0 1 0 0],...  
         'offset',0,'duration',0)
```

Classifier prediction:

```
ev=struct('sample',123,...  
         'type','prediction',...  
         'value',[-1 -1 -1 1 -1],...  
         'offset',0,'duration',0)
```

Imagined Movement event:

```
ev=struct('sample',123,...  
         'type','stimulus.move',...  
         'value','left-hand',...  
         'offset',0,'duration',300)
```

Compact notation:

```
s:123,t:'stimulus.flash',v:[0 0 1 0 0],o:0,d:0
```

(key functions) Event manipulation

evt=mkEvent(*type,value,[sample,offset,duration]*)

- make a buffer event, with sensible defaults

sendEvent

- *evt*=sendEvent(*type,value,[sample,offset,duration,host,port]*)
- *evt*=sendEvent(*evt,[host,port]*)
- Send event to the buffer on machine host at port.

mi=matchEvents(*evts,mtype,mval*)

- Find events with type *mtype* and value *mval* in *evts* a vector of event structures.
- *mtype* – can be cell-array of types to match, e.g. {'type1' 'type2'}
- *mval* – can be cell-array of values to match, e.g. {'val1' 10 'val3'}
- match if any *mtype* matches and any *mval* matches,
 - i.e. above matches (t:'type1',v:10), (t:'type2',v:10),(t:'type1',v:'val1')
- *mi* is logical vector of which *evts* matched
- N.B. Empty ([]) or '*' *mtype/mvalue* matches everything

(key functions) Event manipulation

```
import Fieldtrip
```

```
evt=Fieldtrip.Event(type,value,[sample,offset,duration])
```

- make a buffer event, with sensible defaults

```
sendEvent
```

- *evt=bufhelp.sendEvent(type,value,[sample,offset,duration,host,port])*

```
evtfilter=bufhelp.createeventfilter((mtype,mvalue))
```

- Create a function to filter a list of events and return only events with with type *mtype* and value *mval* in *evts* a vector of event structures.
- *Can also use a list of (mtype,mvalue) pairs to match multiple types, e.g. [('type1','val1'),('type2','val2')]*
- *Or a list of strings to only match on the type: e.g. ['type1','type2']*
- *Apply the resulting filter function:*

```
mevts=evtfilter(evts)
```

Hands-on 1: Event Echo-Server

Experiment Task

- Write a simple echo-server which:
 - Connects to the Buffer server
 - Waits for **any** incoming event, and
 - Responds by sending a 'echo' event with the same value but type='echo'
 - Quits if it receives an event with type='exit'
- N.B. Don't 'echo' your own echo events!

Assignment:

- start from : `echoServer_skel.{m/py}`

N.B. send 'keyboard' events by pressing keys in the eventViewer window.

Useful Functions:

- Setup MATLAB paths and connect to buffer:
 - See the header code block in echoServer_skel.m
 - This will initialize the paths:
 - `run ../../utilities/initPaths.m`
 - Then try to connect to the buffer server until valid header is returned.
 - `while (isempty(hdr))|.....`
 - Then initialize some utility functions for high-precision timing
 - `initsleepSec();initgetTime();`

Useful Functions:

[devents,state]=...

buffer_newevents(*host,port,state,mtype,mval,timeout_ms*)

- **wait** for any **new** events **matching** (*mtype,mval*)
 - Matching done by *matchEvents*
 - *mtype* – can be cell-array of types to match, e.g. {'type1' 'type2'}
 - *mval* – can be cell-array of values to match, e.g. {'val1' 10 'val3'}
 - match if **any** *mtype* matches **and any** *mval* matches
- return the matched events in the vector of structure(s) *devents*
- *state* is the match state, used to track which events have been processed between function calls
- Return after *timeout_ms* milliseconds even if no matching events found

Useful Functions:

- Setup PYTHON paths and connect to buffer:
 - See the header code block in echoServer_skel.py
 - This will initialize the paths:
 - `bufhelpPath = "../..python/signalProc"`
 - `sys.path.append(os.path.join(os.path.dirname(os.path.abspath(__file__)), bufhelpPath))`
- Then try to connect to the buffer server until valid header is returned.
 - `ftc, hdr = bufhelp.connect();`

Useful Functions:

`devents=bufhelp.buffer_newevents(mtype,timeout_ms)`

- **wait** for any **new** events **matching** (*mtype*)
 - Matching done by *matchEvents*
 - *mtype* – can be list of type strings to match, e.g. ['type1', 'type2']
 - match if **any** *mtype* matches **and** **any** *mval* matches
- return the matched events in the vector of structure(s) *devents*
- *state* is the match state, used to track which events have been processed between function calls
- Return after *timeout_ms* milliseconds even if no matching events found

Debug/Test your echoServer

- **Start an event viewer:** `dataAcq/startJavaEventViewer.{sh,bat}`
- Or directly buffer+simulated EEG+eventViewer: `debug_quickstart.{sh,bat}`
- **Generate events to be echo'ed:**
 - Directly in the javaEventViewer.
 - Press return
 - Enter the event type, e.g. 'hello'
 - Press return
 - Enter the event value, e.g. 'world'
 - Press return
 - You should now see a copy of your created event in the event log.
 - If the your echoServer is working you should also see your 'echoed' version.

Echo-Server in different languages

Basic echo-server example has been implemented in multiple languages

- **MATLAB:** `echoClient/matlabclient.m`
- **JAVA:** `java/echoClient/javaclient.java`
- **C#:** `csharp/echoClient/csharpclient.cs`
- **Python:**
`python/echoClient/pythonclient.py`
- **C:** `c/echoClient/cclient.c`

Hands-on 2: “Hello World”

Experiment Task

- Display the string “Hello World” (or any other pre-specified string) on the screen, and wait for a key to be pressed to exit
- Send events to annotate what has happened, e.g. startup, string display, key-pressed, shutdown etc.

Method:

- Start from the 'helloworld_skel.{m/py}' function skeleton
 - contains initialisation code to connect to the ft_buffer
 - Some examples of functions you may find useful

Note : event timestamps

- Accurate event time-stamps are **critical** for **evoked** potential analysis
 - >10ms event jitter causes significant reduction in signal quality
- However,
 - data-acquisition may only send data every >20ms
 - And this data may be subject to additional network delays of >20ms
- Stop this **jitter** reducing time-stamp accuracy by;
 - ~~aligning (and tracking) computers real-time-clock and data-sample clock to prevent this jitter reducing~~

Hands-on 3: Sequenced Sentences

Experiment Task

- display set of sentences on the screen where every second 1 more character gets added to the sentence
- pause for 5 seconds between sentences (and/or wait for key press)
- send events for everything that happens

Assignment:


- Make flowchart
- Write code -> test -> debug -> until it works :-)
- Start from `runSentences_skel.{m/py}`

Useful Functions:

- MATLAB: `sleepSec(time)` sleep for the indicated duration in second
- PYTHON: `sleep(time)` sleep for the indicated duration in seconds

Flowchart : sequenced sentences


sequenced sentence

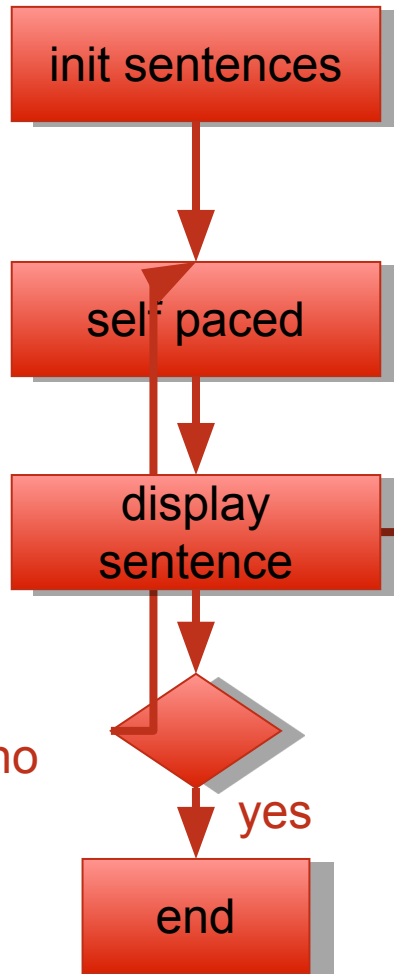

next
sentence

end of sequence

Sentence

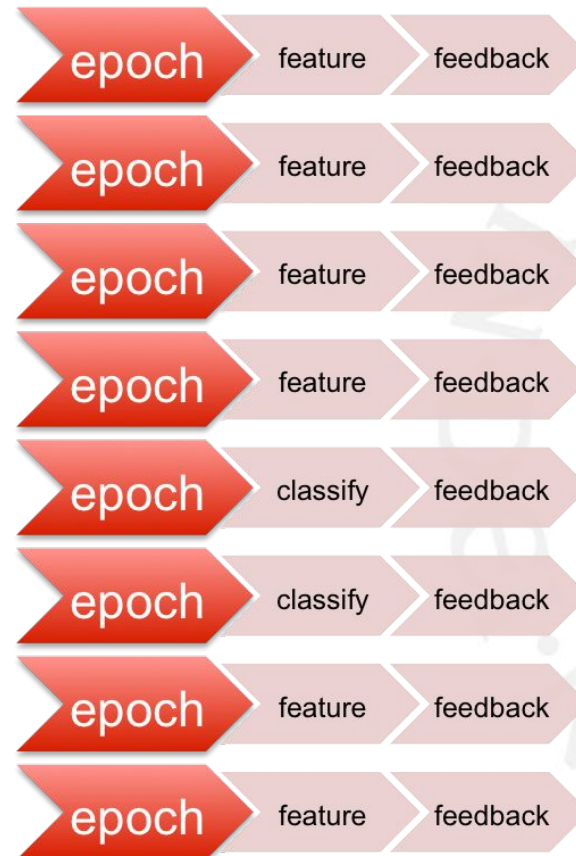


Stimulus presentation


M
a
k
i
n
g
B

.Events and processing functions

sequenced sentences

Events

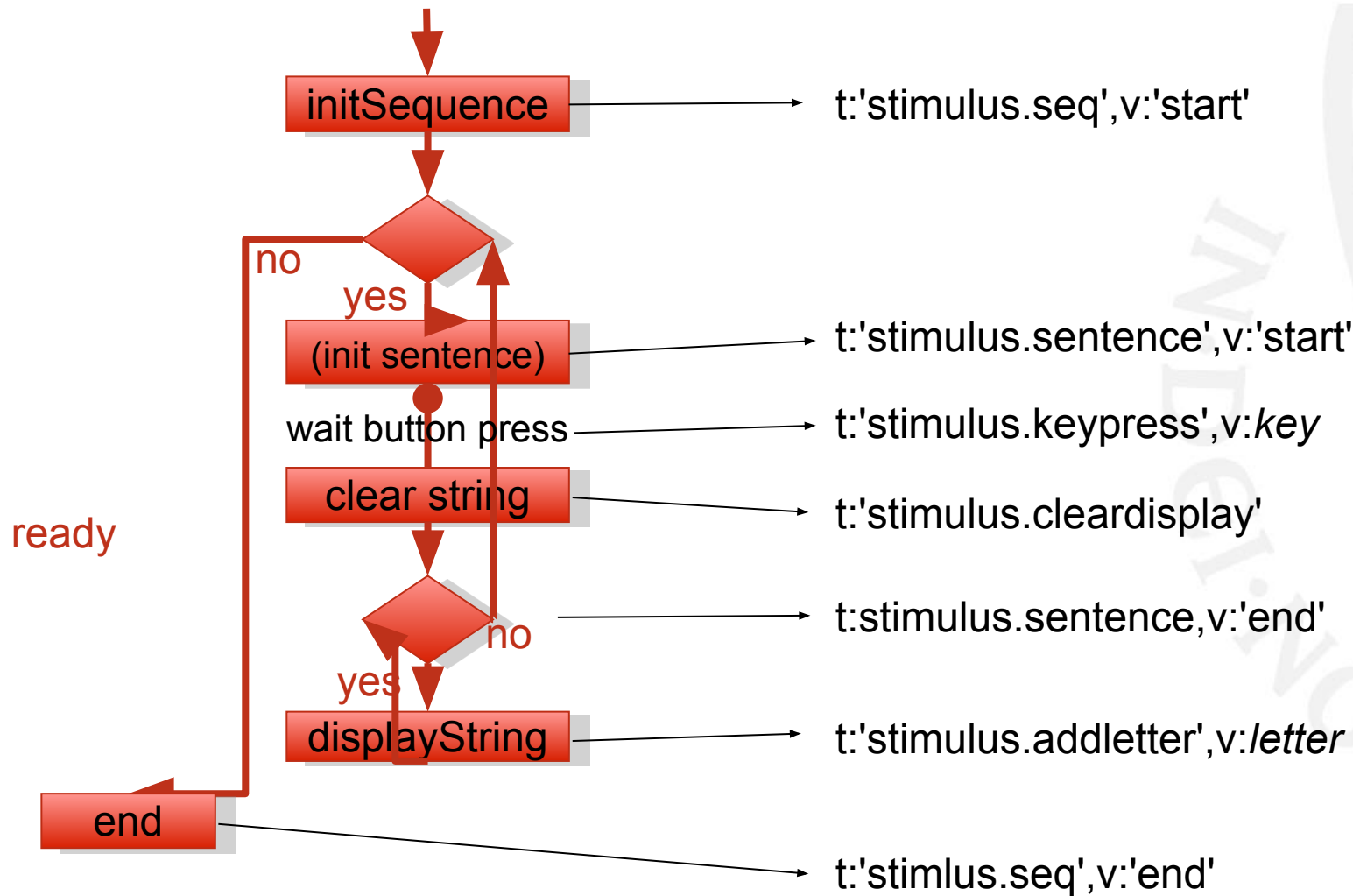

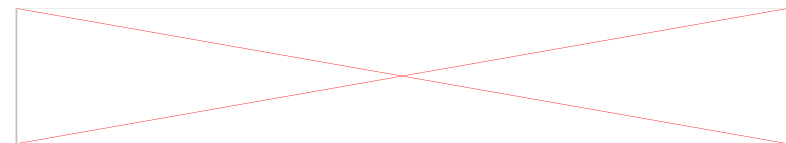

Hands-on 4: ERP Viewer

Experiment Task

- In 5 sequences of 10 seconds:
 - Every 1 seconds: either randomly display or don't display a cross (+) on the screen for 200ms
- Display a 'Press key to continue string' between sequences, and wait for key press to move to the next sequence
- For every 'stimulus event', i.e. point when the '+' could have been displayed, record 600ms of data annotated with whether it was a '+' or not
- Every time you get some data, compute an average of the EEG data for that type of stimulus, i.e. + or no-+, and display the resulting averages as a multi-plot on the screen
-
- N.B. You will need a separate signal processing process to: get the data, compute the ERP and display the results!
- Assignment:
- Make flowchart for each of the processes, i.e. stimulus, and signalProcessing
- For the expt-control & stimulus presentation start from : `runStimulus_skel.{m/py}`
- For the signalProcessing & results generation use : `runSigProc_skel.{m/py}`

Useful Functions:

[data,devents,state]=buffer_waitData(*host,port,state,...*

'startSet',startEvts,'trlen_samp',samp,'exitSet',exitEvts)

- **for all** events **matching** *startEvts* record *samp* samples of data
- **until** an event matching *exitEvts* is generated
- *startEvts* and *exitEvts* specify the events to match in the format:
 - type – event type has this value
 - {'type' val} – event has type==type and value==val
 - {'type1' 'type2'} – event has type == 'type1' or 'type2'
 - {'type1' 'type2'} {val1 val2} – event has type == 'type1' or 'type2' **and** value==val1 or val2
- return the matched event structure(s) in *devents* and corresponding data in *data*
 - *Data* is a vector of structures. data.buf = [nChannels x nSamples] raw EEG data
- *state* is the match state, used to identify which events have been processed between function calls
-
- N.B. *ExitEvts* has the special event type 'data' which returns as soon as the data is available for the first matched *startEvt*

Useful Functions:

`[data,devents,stopevents]=bufhelp.gatherdata(trigger,trlen,exitTrigger)`

- **for all** events **matching** *trigger* record *trlen* samples of data
- **until** an event matching *exitTrigger* is received
- *trigger* and *exitTrigger* specify the events to match in the format used for `bufhelp.createeventfilter`. e.g.
 - `[type1, type2,...]` - list of types to match
 - `[(type1,val1),(type2,val2)...]` - list of type+val pairs to match
- return the matched event structure(s) in *devents* and corresponding data in *data*
- *Data* is a vector of structures. `data.buf = [nChannels x nSamples]` raw EEG data

DEBUGGING: ERP Injection

- Debugging the correct signal-analysis script is difficult without a **true-signal** to validate that the trigger event timing/value is correct.
- To make this easier the **simulated-eeg** data supports **ERP-injection**.
 - This allows the presentation system to ‘tell’ the EEG simulator to add an erp to the **TRG** channel at this point.
 - This can then be used to check the correctness of the signal processing, e.g. do all the ‘flash’ events have a trigger ERP at time=0? do all the ‘non-flash’ events have no ERP?

DEBUGGING: ERP Injection

- ERP injection works by sending a number for the ERP strength on UDP port 8300.
- Setup the socket connection:

```
trigsocket=javaObject('java.net.DatagramSocket');  
trigsocket.connect(javaObject('java.net.InetSocketAddress','localhost',8300));
```
- Use the socket to inject an ERP

```
trigsocket.send(javaObject('java.net.DatagramPacket',int8([1  
0]),1));
```

DEBUGGING: ERP Injection

- ERP injection works by sending a number for the ERP strength on UDP port 8300.

- Send a number to inject an ERP:

```
import socket
```

```
socket.socket(socket.AF_INET,socket.SOCK_DGRAM,0).sendto(bytes(1),('localhost',8300))
```

Summary

- BCI can be broken into 4 processes: data-acquisition, experimental control, signal processing, and stimulus presentation
- `buffer_bci` framework : uses buffer events as a blackboard for inter-process communication